

Smart Chart system with Full HD LCD
PLC-8000, PLC-8000pola



SPECIFICATION

Test Distance	1.0~6.0m (Unit:0.1m)
Chart	Randot Ring, Alphabet, Numerals, Snellen, Pictures, Mask (Vertical, Horizontal, Single letter), RG chart/ Polaroid chart (PLC-8000pola)), Color blind, Amsler Chart
Programmable chart	Up to 40 charts
Interface	USB(USB2.0),RS-232
Communication	Wired : RS-232 Wireless : Bluetooth
LCD	23 Inch Full-HD (Resolution : 1920x1080),PLC-8000 24 Inch Full-HD (Resolution : 1920x1080),PLC-8000pola
OS	Android 4.4.2
Storage	Micro SD Card(8G)
Speaker unit	1.5W
Video format	AVI, MPEG, MP4, MKV
Rated Voltage	AC100-240V, 50/60Hz
Power consumption	PLC-8000 : 45-55VA PLC-8000pola : 35-45VA
Volume (Dimension)	Body : 583(W)×55(D)×423(H) mm Remote controller : 64(W)×196(D)×20(H) mm
Weight	Body : 약 6kg, Wall MountBracket : 370g Remote control : 160g (Battery included)
Environment conditions	(In use) Temperature : +10 ~ +40 Humidity : 30%RH ~ 85%RH Pressure : 700hPA ~ 1060hPA
	(In transport/storage) Temperature : -10 ~ +55 Humidity : 10%RH ~ 90%RH Pressure : 600hPA

PT-0016613

⁶Designs and details can be changed without prior notice when its improvement.

POTEC POTEC Co., Ltd.
40-4, Techno 2-ro, Yuseong-gu, Daejeon 34015, Korea
Tel +82 42 632 3536 Fax +82 42 632 3537
webmaster@potec.biz www.potec.biz

Distributde by
BIRMINGHAM OPTICAL
Unit 4 Gravelly Industrial Park, Birmingham B24 8HZ
T: 0808 123 2020
E: sales@birminghamoptical.co.uk | www.birminghamoptical.co.uk



POTEC

Make smart technology yours



PLC-8000
PLC-8000pola

Smart Chart system with Full HD LCD




PLC-8000 PLC-8000pola

Smart Chart system with Full HD LCD

CONVENIENCE

- Smart Chart system with Full HD LCD
 - 23(Regular), 24 Inch (Polaroid) High resolution chart
 - Full-HD TFT-LCD(1920x1080) panel / 16,7MColor 250cd/m²
- Wide range of test
High Resolution LCD enables wide range of test 0.05 ~ 2.0
- Adjustable test distance
Test distance is adjustable 1.0m ~ 6.0m by 0.1m
- SMART intercommunication
Support both interface RS-232(Wired), Bluetooth(Wireless) for the compatibility with PDR-7000
- Basic test / User modified program
Provide basic Program A and Program B which provides programmable up to 40 steps.
- Enhanced R/G filter adjustment
Broaden range of R/G adjustment up to 96 steps enables optometrist to optimize the screen to fit in well with environment for reliable binocular vision test.
- Promotional Videos Play
During the power saving mode, the LCD backlight is turned off to minimize power consumption and the promotional videos for clinics and optical shops are playable on a LCD chart.
- User friendly functions
Contrast, Mirror, Mask.
- Customizable boot-up screen & slides
(Up to 38 slides)

OPTOMETRY

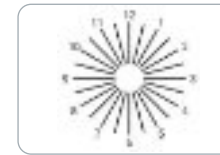


Diversified Astigmatism test chart
Expanded diagnose environments with various charts (7 types)



Astigmatism Fan & Block Charts
The addition of fan and block charts to general astigmatism test functions allows more accurate astigmatism test.



Reverse mode for Lasik and Cataract surgical patients

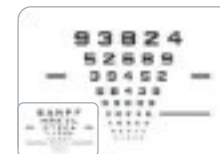

Diversified ETDRS test chart
Newly added letter & number based charts provide wide range of diagnose environments.



Amsler Vision test chart for Macular Degeneration

Enhanced Maddox rod test



Embedded Fixation Target Lamp for Maddox test provides 10 level of brightness

Full-HD Video play



Playing cartoon film at High resolution helps you draw attention from kids or infants that are easily distracted for intensive diagnose. Play full HD video from USB flash drive and control by remote controller.

Various color blindness test chart

